


Handicare Simplicity

A Minimum swivel radius – with short arm	680mm 644mm
B Minimum wall to seat back	75mm
C Seat depth	370mm
D Minimum folded width to footplate	285mm
E Minimum open width to edge of footplate	580mm
F Armrest width – external	625mm
G Armrest width – internal	480mm
H Seat back height	400mm
I Footplate to seat height – standard seat option – extended seat post	390-490mm 460-560mm
J Minimum footplate height	50mm
Bottom parking space required	660mm [†]
Minimum track intrusion into staircase	160mm
Minimum track intrusion into top landing (typical at 45°)	150mm*
Maximum weight limit – Hinge/Slide Track	140kg/22 stone 120kg/19 stone
Speed	0.12m/sec (approx)
Power supply	24v DC battery
Charger	240v
Staircase angle**	27-50°

[†]Increases by up to 140mm with zero intrusion kit.

*A 0mm track intrusion kit is available.

**Recommended minimum Slide Track staircase angle 37°.

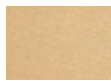

BUTE



handicare

Handicare Simplicity

Straight stairlift



The benefits of a Handicare Simplicity straight stairlift

If you are thinking about installing a stairlift the Handicare Simplicity offers you a safe and reliable solution.

Handicare stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

The Handicare Simplicity straight stairlift offers you all the basic functions that you would expect from a stairlift and, if you have a doorway close to the base of the stairs, you can choose from a manually operated hinged track or Handicare's unique

Slide Track solution to ensure the doorway can be left clear of obstructions when your stairlift is not in use.

Supplied in a neutral shade, the wipe clean upholstery will compliment most colour schemes.

The track is fitted to the stairs, not the wall, so it is unlikely you will need to redecorate following the installation of your Handicare Simplicity stairlift.

Like all Handicare stairlifts, the Simplicity has a battery backup which means it will continue to work in the event of a power cut so you do not need to worry about being stranded up or downstairs.

The compact unit folds up leaving the staircase free for other users and the Handicare Simplicity also comes with two remote controls which means, once folded, you can move it along the track freeing up space in your hallway.



Safety as standard

Your Handicare Simplicity stairlift will come with a range of safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs.

Your Handicare Simplicity stairlift has a clip operated seat belt for added security and peace of mind.

Easy to use

The operating toggle is easy to use, simply hold it towards the direction you wish to travel, up or down, and your Handicare Simplicity will move smoothly and gently on its way.

The control is conveniently located at the end of the arm leaving both arms clear to provide unhindered support when getting on or off the seat.

